



VEBKO ONLINE HELP

“New Based on Custom Template” Option

Category	tester
Section	Modules Toolbar and Statusbar
Language	English
Generated	2026-10-18

In addition to "New" option, this option is designed to ease the test and save time. It is possible for you to create a default "Template" in accordance with your needs and when necessary, by loading this "Template" apply the saved changes.

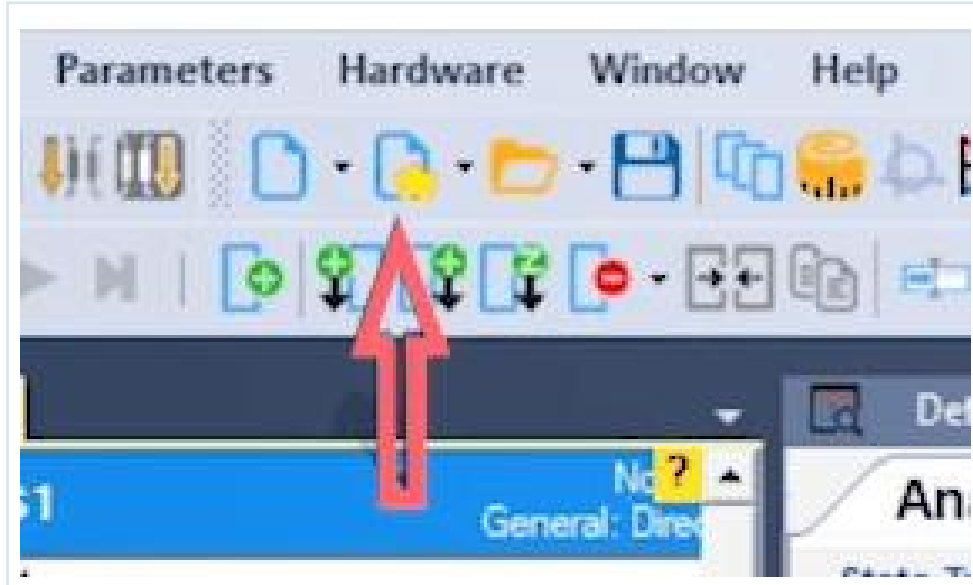


Figure 1

1 First Mode: "Custom Template"

In "Custom Template" you can create a "Template" by using "Device Template", "Report Template", or "Decimal Places Template" modes. This can be explained by creating a "Template". To do so, in "Preferences" in "Room" section, you can create a "Template" from "Device Template" and "Template Report Package" or "Decimal Places Template".

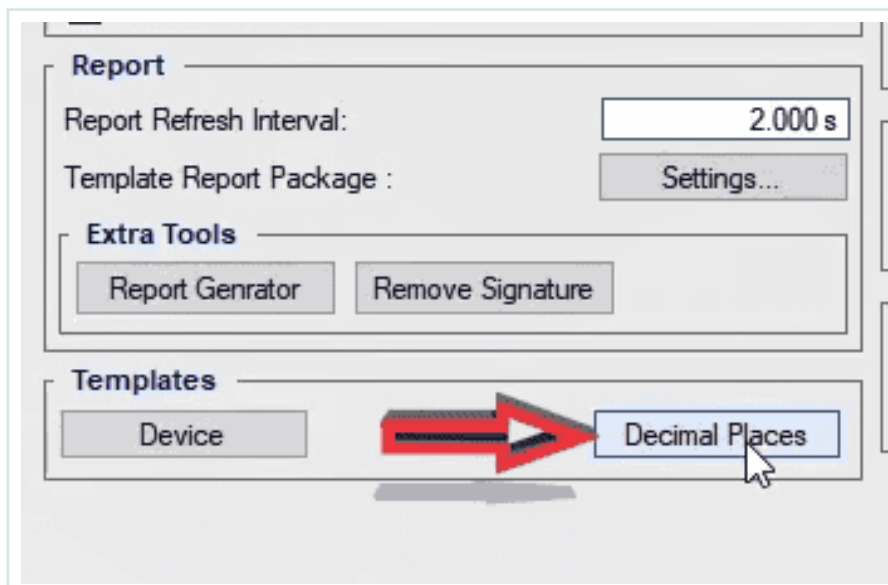


Figure 2

For example, after entering "Preferences" window, in "Decimal Places Template" section for "Sequencer" room, select and apply "Volt" as the unit for voltage in "Template". Enter "Vebko" in "Name" and "Manufacturer" fields in "Device" section, and then by applying the settings, save these changes as a "Template". Now, to use these "Templates" in the rooms, you need to select the saved "Templates". If you open "Decimal Places" in "Sequencer" room, you can see that the changes made in "Preferences" is not applied to this section. Now to apply the "Template", by clicking on "New Based On Template" in "New" window, you can see "Custom Template" and "Saved Template" radio buttons. By selecting "Custom Template" radio button, you can specify which of the three mentioned sections is to be applied to the new file. Then, by confirming the settings by selecting "Ok" you can view the applied changes.

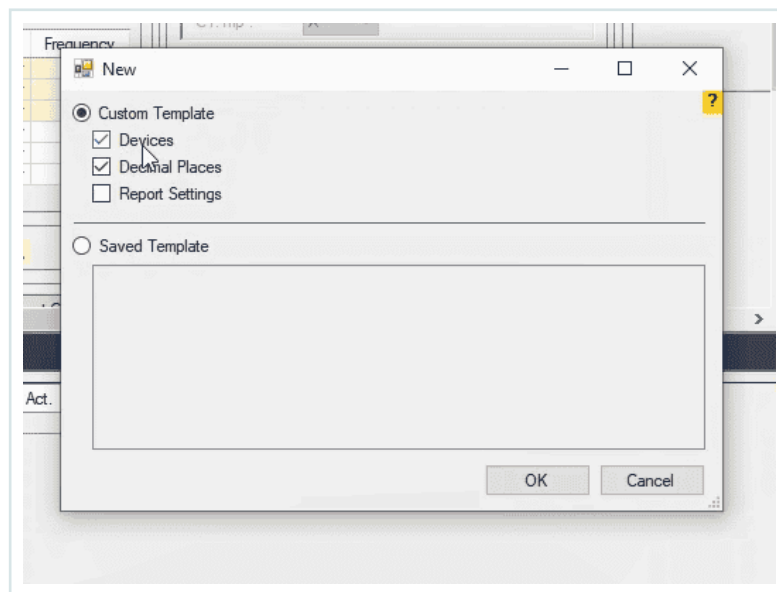


Figure 3

2 Second Mode: "Saved Template"

In "Saved Template", you can specify the default necessary windows, arrangements and the number of fields. To better understand this, select "State Type Continuous", open "Vector View" window and specify the desired size. Then, select "Save as Template" from 'File' menu and select a title for it. If you wish to apply your specified "Template" in a new room, you need to click on "New Based On Template" and select "Save Template" radio button. In the end, select the desired "Template" and apply it.

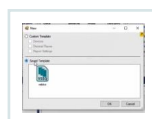


Figure 4

[jg1]New_Based_On_Custom_Template_Option